



Contact Gown



N95 Respirator



Face shield or goggles



Gloves



## Gown + Gloves

1 **Remove** any personal items and jewelry and put in secure location, not in pockets.

2 **Sanitize** hands.

3 **Put on** contact gown outside room.

Open-end **faces** your back.

**Tie** the back of the gown.

4 **Put on** gloves over the cuffs of the gown.

## Respirator

5 **Put on** N95, ensuring proper seal. **Ensure** straps are not crossing.

6 **Place** hands over the front of the N95. **Breathe** an easy deep breath in and out.

If you **feel** air escape the edges, **refit** and **repeat**.

7 **Sanitize** gloves.

## Eyes + Entry

8 **Put on** face shield.

9 **Sanitize** gloves.

10 **ENTER** room

**DO NOT** enter the room if you do not achieve a proper respirator seal.

11 **Do not** touch face or re-adjust N95 or face shield inside room.



### ACE, contact gown, gloves

## Gown + Gloves


1 **Sanitize gloves.**

2 Cross arms and **grip** gown on shoulders. **Pull and break** gown in controlled fashion.

3 **Roll** the gown towards your hands. **Remove** the gloves with the gown. **Dispose** of gloves and gown.

4 **Sanitize** hands.

**EXIT patient room**

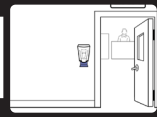

**EXIT patient room**

## Eyes

5 **Put on** new gloves.

6 **Sanitize** gloves.

7 **Do not** touch face.

8 **Remove** face shield by the **strap** over your head without touching your skin.

## Mask or Respirator

9 **Sanitize** gloves.

10 **Pinch** bottom strap and **pull** far over head. **Do not** let straps touch your face.

11 **Pinch** top strap and **pull** far over head. **Do not** let straps touch your face as you **remove** the N95.

## Wash

12 **Remove** gloves.

13 **Head immediately** to handwashing station. **Wash** hands with soap and water.